

Introduction

Anchored Multiplex PCR (AMP™) assays, VARIANTPlex and FUSIONPlex are amplicon-based, next generation sequencing (NGS) workflows that identify genomic variations present in DNA or RNA input, respectively. AMP reagents are available in a lyophilized format which are user-friendly and convenient to store. However, lyophilized reagents are not suitable for automated liquid handling (ALH) or high-throughput (-HT) environments. Thus, we developed AMP reagents in liquid format to satisfy the need for a platform-agnostic method of high-throughput AMP library preparation. Using an ALH platform we developed, tested and optimized workflows for VARIANTPlex-HT and FUSIONPlex-HT liquid reagents. Taken together, we provide a solution for increasing throughput for AMP-based NGS customers while reducing hands-on time.

Methods

The performance of our liquid reagent based VARIANTPlex-HT and FUSIONPlex-HT workflows were evaluated relative to legacy lyophilized reagents using several inputs of varying quantity, quality and type (RNA & DNA). We examined liquid reagent performance across multiple AMP panels, which target solid tumor and blood cancers variants. We also demonstrate the relative performance of the automated FUSIONPlex-HT and VARIANTPlex-HT workflows to manual library preparation using SeraSeq® FFPE and Myeloid input materials with AMP-based panels. Samples were sequenced on the Illumina NextSeq® or MiSeq® system and analyzed using Archer™ Analysis v7.0 software. Two-sided, unpaired Student's t-Test was used to indicate a statistically significant difference between means. For multi-group comparisons against a single control condition, the Dunnett test was performed. *P<0.05; **P<0.01; ***P<0.001; n.s., not significant.

New! Archer™ Liquid Reagents



- Available in 24- and 96- reaction kits
- Reagents are color matched to protocol steps for ease of use
- Each reagent supplied with 20% overage
- Compatible with all VARIANTPlex and FUSIONPlex panels
- Features new liquid adapters
 - 96-well plate format
 - 192 unique P5 MBC + P7 adapters for a total of 36,864 unique combinations!
 - Reduced PhiX Control requirement to 1-2%

Results - FUSIONPlex-HT

Input	Input Mass	Panel Name	Panel Size	Depth
Seraseq Fusion RNA Mix v4	20ng	FUSIONPlex Lung v2	17 genes	1M
Seraseq Fusion RNA Mix v4	50ng	FUSIONPlex Lung v2	17 genes	1M
Seraseq Fusion RNA FFPE v4	50ng	FUSIONPlex Pan ST v2	137 genes	3.5M

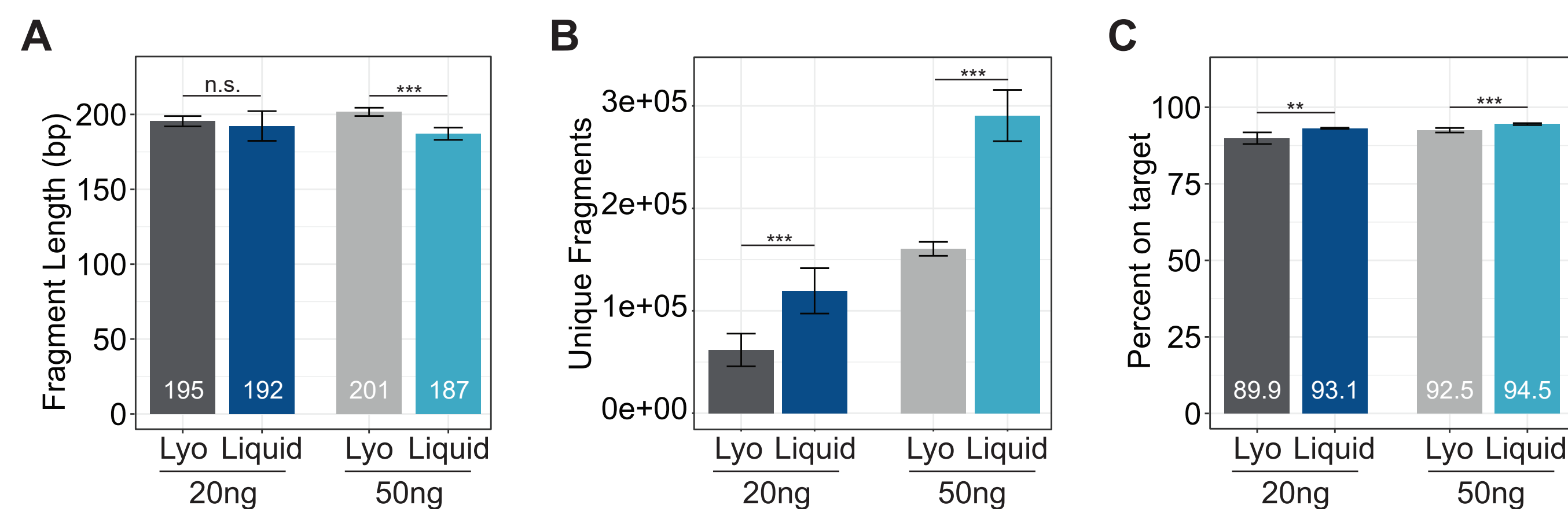


Figure 1. FUSIONPlex Lung v2 panel performance with 20 or 50ng Serare Fusion RNA v4 input. A) Mean fragment length. B) Mean number of unique reads passing read cleaning. C) Mean percent of unique reads that mapped on target. Data shown is aggregated between 2 users (n=3 technical replicates per condition per user). All samples were subsampled to 1M reads.

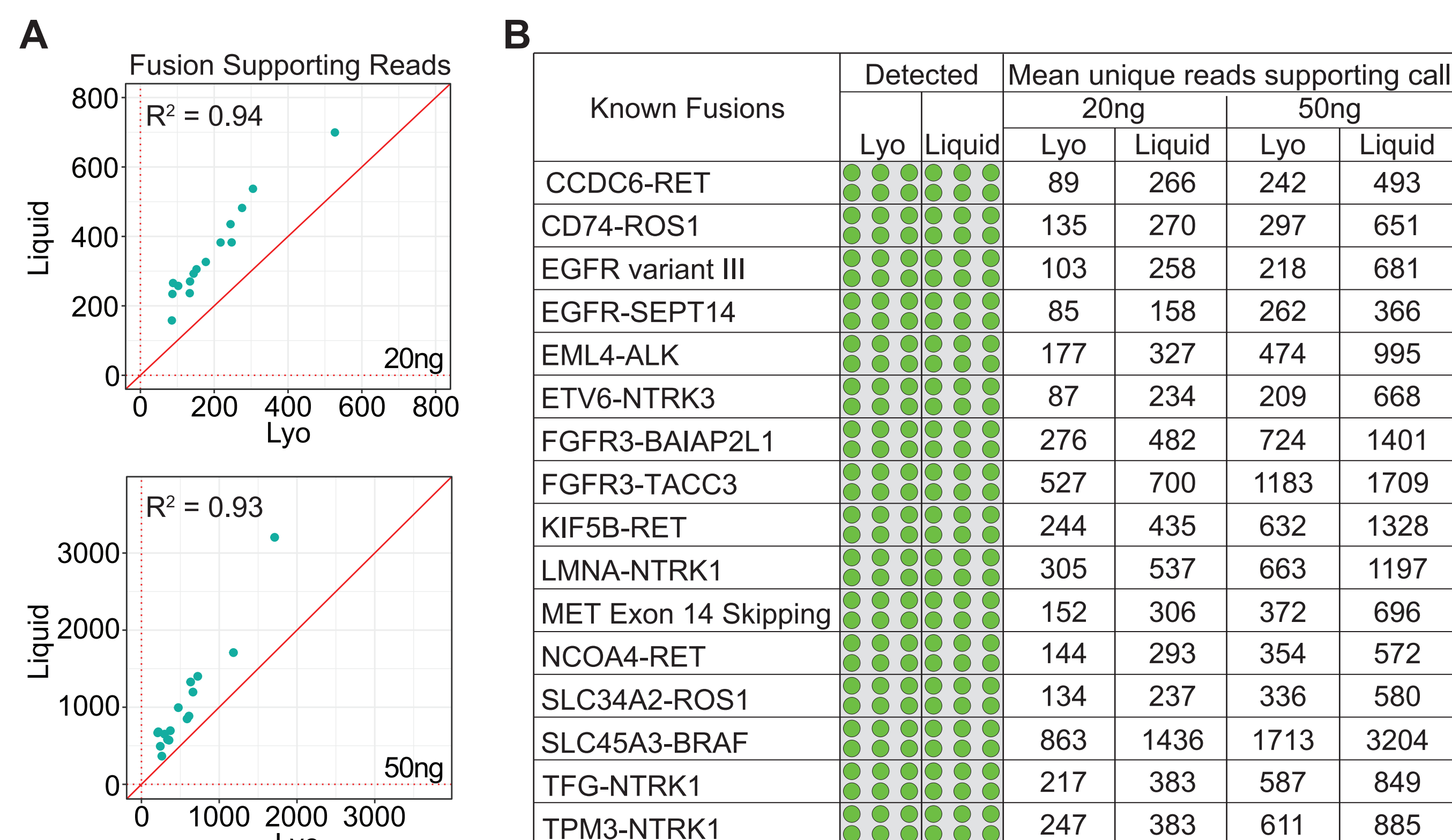


Figure 2. Unique reads supporting known fusions targeted by FUSIONPlex Lung v2 panel with 20ng or 50ng Serare Fusion RNA v4 input. A) Correlation of unique fusion supporting reads between Liquid (y-axis) and Lyo (x-axis) workflows for each known fusion (n=11) using 20ng (Top) or 50ng (Bottom) B) Mean unique fusion supporting reads supporting call of expected fusions. Each green dot indicates the fusion was detected in each replicate (n=3 technical replicates per condition for each user). All samples were subsampled to 1M reads.

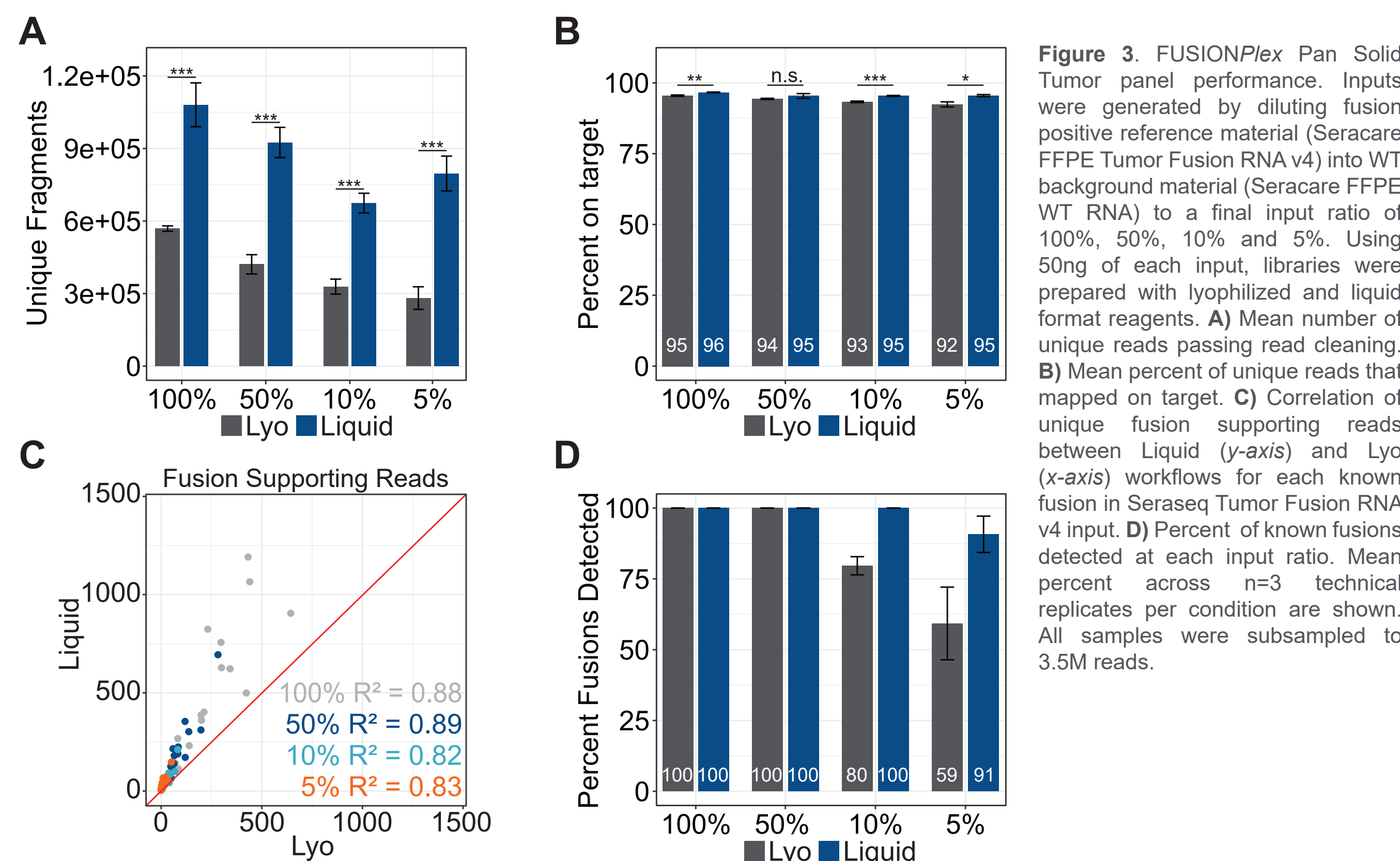


Figure 3. FUSIONPlex Pan Solid Tumor panel performance. Inputs were generated by diluting fusion positive reference material (Serare FFPE Tumor Fusion RNA v4) into WT background material (Serare FFPE WT RNA) to a final input ratio of 100%, 50%, 10% and 5%. Using 50ng of each input, libraries were prepared with lyophilized and liquid format reagents. A) Mean number of unique reads that mapped on target. B) Mean percent of unique reads that mapped on target. C) Correlation of observed allele frequency (AF) between Liquid (y-axis) and Lyo (x-axis) workflows for each known fusion in Serareq Tumor Fusion RNA v4 input. D) Percent of known fusions detected at each input ratio. Mean percent across n=3 technical replicates per condition are shown. All samples were subsampled to 3.5M reads.

Results - VARIANTPlex-HT

Input	Input Mass	Panel Name	Panel Size	Depth
Seraseq Myeloid DNA	50ng	VARIANTPlex Core Myeloid	17 genes	3M
Seraseq Compromised FFPE DNA	50ng	VARIANTPlex Complete ST	137 genes	50M
Horizon Moderate FFPE DNA	10ng	VARIANTPlex Core ST	60 genes	4.5M
Horizon Moderate FFPE DNA	50ng	VARIANTPlex Core ST	60 genes	4.5M

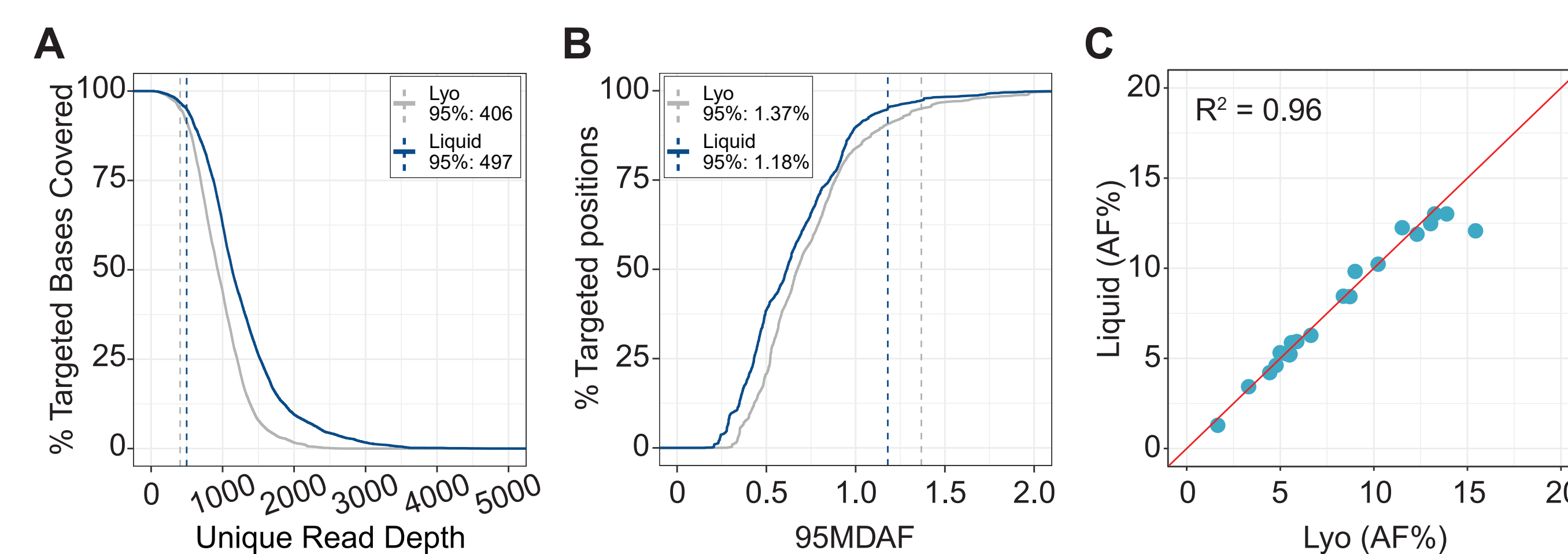


Figure 4. VARIANTPlex Core Myeloid panel performance with 50ng Serare Myeloid DNA mix input. A) Mean unique read depth at bases targeted by the panel. Upper right: mean unique read depth 95 percentile. B) Minimum detectable allele frequency (MDF) at alpha = 0.05 across bases targeted by the panel. Upper left: 95MDF 95 percentile. C) Correlation of observed allele frequency (AF). 151/154 (97%) and 154/154 (100%) of expected variants were observed across Lyo and Liquid workflows, respectively. All data shown is aggregated between 2 users with at least n=3 technical replicates per condition per user. All samples were subsampled to 3M reads.

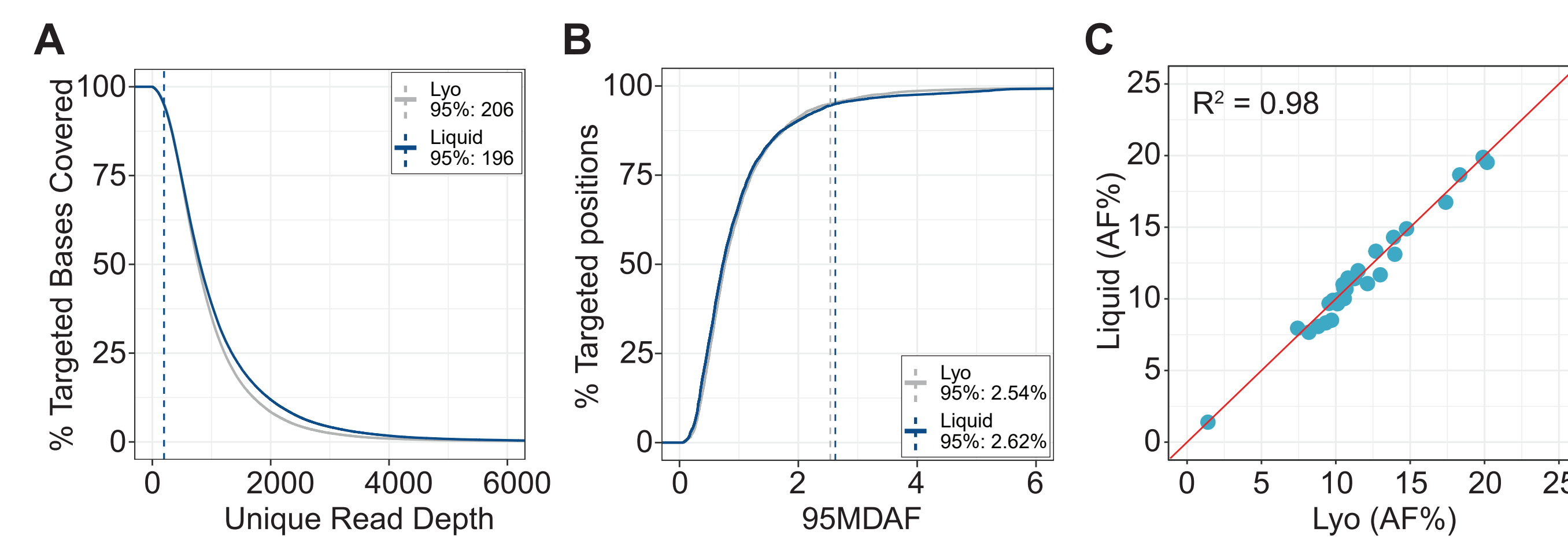


Figure 5. VARIANTPlex Complete Solid Tumor panel performance with 50ng Serare compromised FFPE DNA input. A) Mean unique read depth at bases targeted by the panel. Upper right: mean unique read depth 95 percentile. B) Minimum detectable allele frequency (MDF) at alpha = 0.05 across bases targeted by the panel. Lower right: 95MDF 95 percentile. C) Correlation of observed allele frequency (AF). 84/84 (100%) of expected variants were observed across both Lyo and Liquid workflows. All samples were subsampled to 50M reads (n=3 technical replicates per condition).

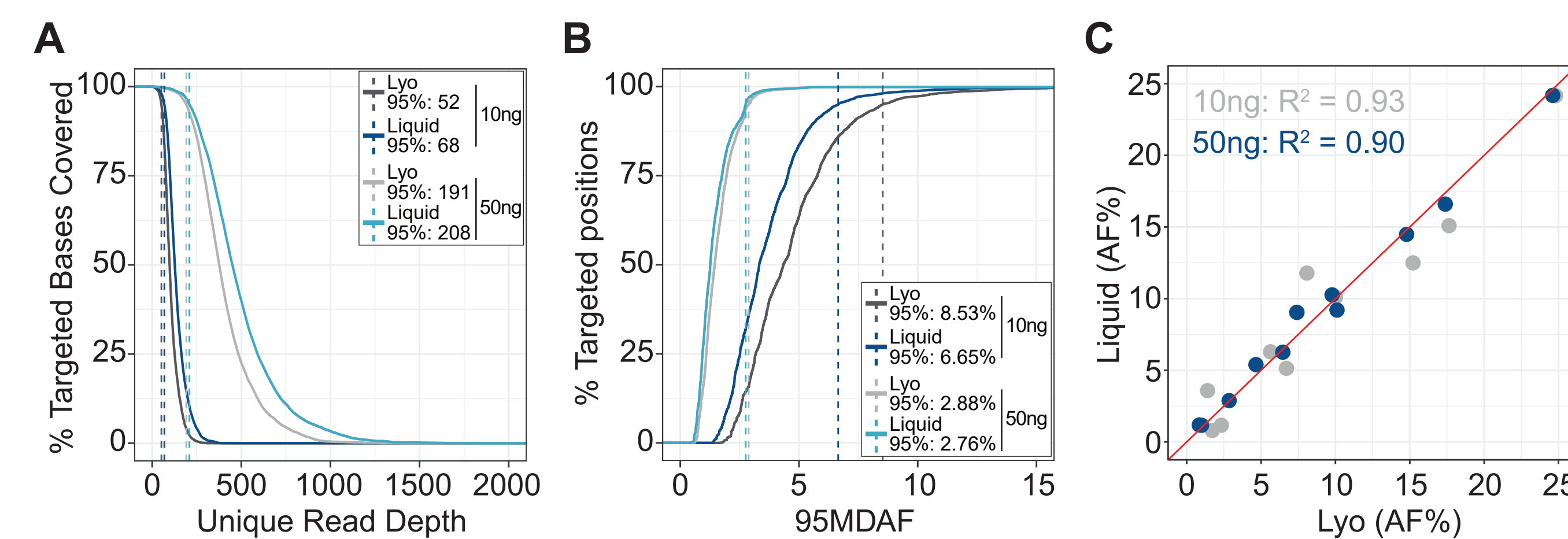


Figure 6. VARIANTPlex Core Solid Tumor panel performance with 10 or 50ng Serare Horizon Moderate FFPE DNA input. A) Mean unique read depth at bases targeted by the panel. Upper right: mean unique read depth 95 percentile. B) Minimum detectable allele frequency (MDF) at alpha = 0.05 across bases targeted by the panel. Lower right: 95MDF 95 percentile. C) Correlation of observed allele frequency (AF). 152/154 (98%) and 163/165 (98%) of expected variants were observed across both Lyo and Liquid workflows. All data shown is aggregated between 2 users with at least n=3 technical replicates per condition per user. All samples were subsampled to 4.5M reads.

Results - Liquid Automation

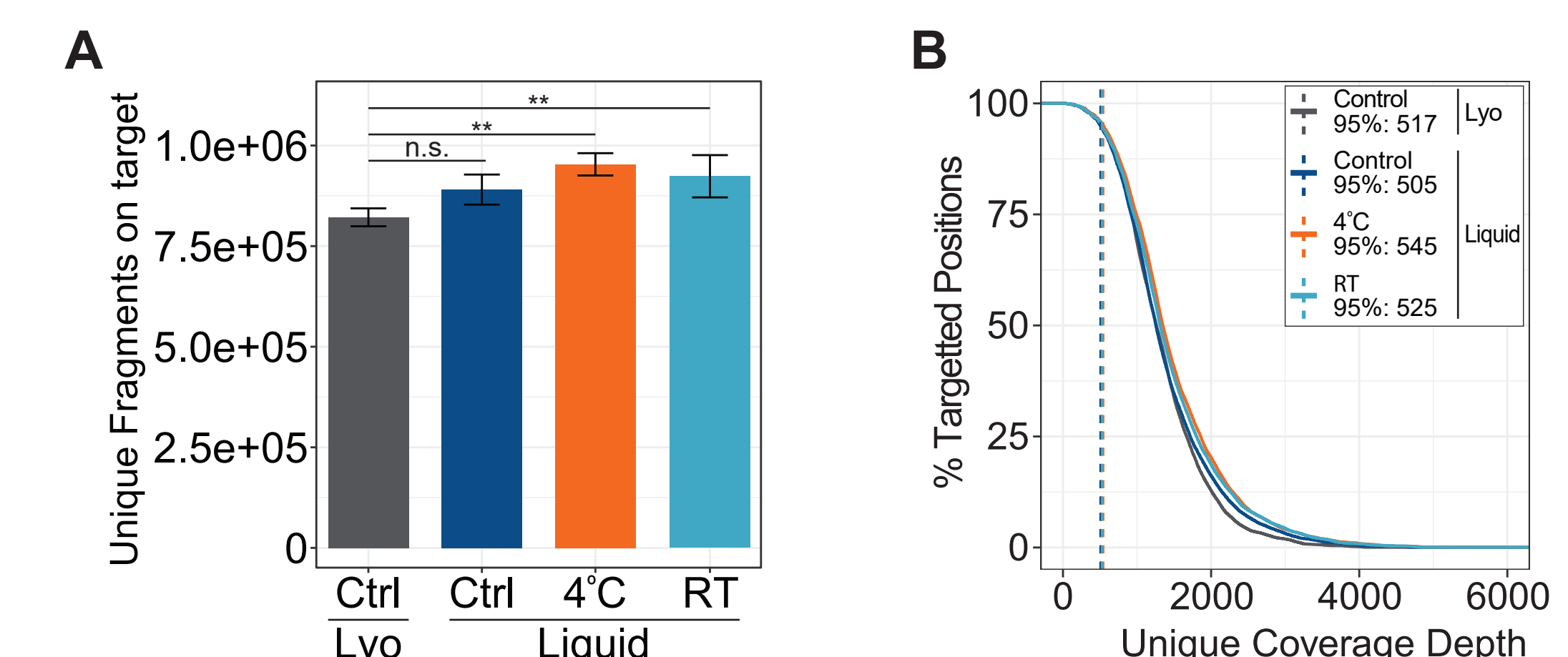


Figure 7. VARIANTPlex Core Myeloid panel performance with 50ng male Genome in a Bottle input. Liquid mastermixes were prepared at the step where required (Liquid Ctrl; Control) or all mastermixes were prepared at the start of the library prep and stored at the designated temperature (4°C or RT; room-temperature) until required at the designated step of the workflow. A) Mean unique reads mapped on target. B) Mean unique read depth at bases covered by the panel. All samples were subsampled to 3M reads (n=3-4 technical replicates per condition).

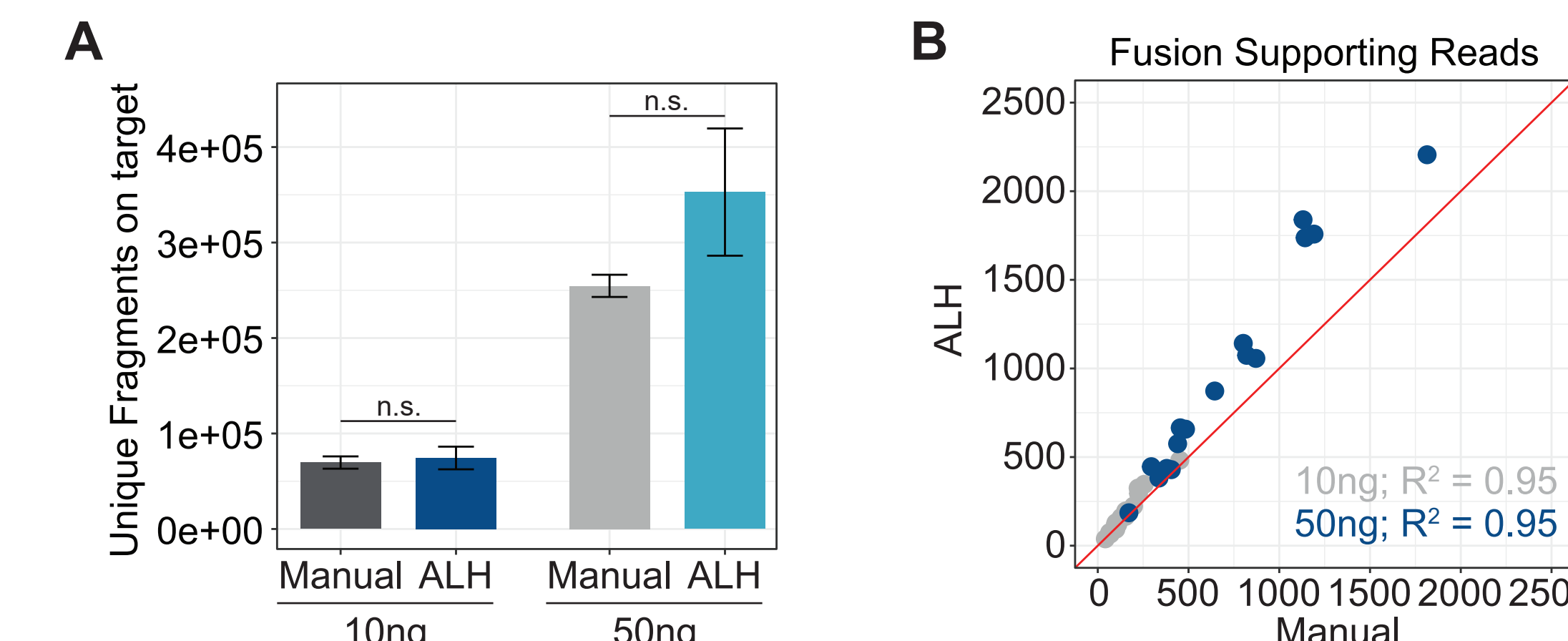


Figure 8. FUSIONPlex Lung v2 panel performance with 10ng or 50ng Serare FFPE Tumor Fusion RNA v4 input. Liquid libraries were prepared simultaneously either manually or on an Automated Liquid Handler (ALH) prior to sequencing. A) Mean unique reads mapped on target. B) Correlation of unique fusion supporting reads between Manual and ALH-prepared samples. 48/48 (100%) of expected fusions were detected by both workflows. All samples were subsampled to 1M reads (n=3 technical replicates per condition). Error bars are mean ± SEM.

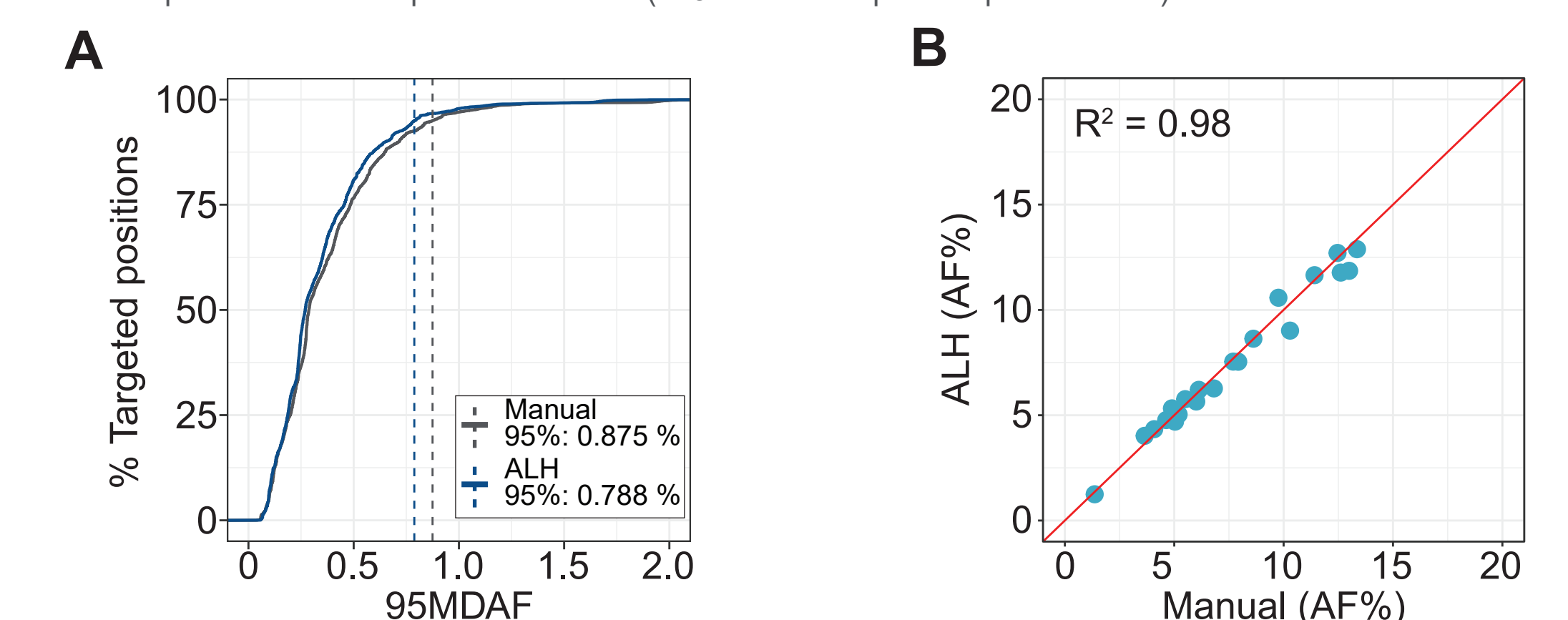


Figure 9. VARIANTPlex Core Myeloid panel performance with 200ng Serare Myeloid DNA mix input. Liquid libraries were prepared simultaneously either manually or on an Automated Liquid Handler (ALH) prior to sequencing. A) Minimum detectable allele frequency (MDF) at alpha = 0.05 across bases targeted by the panel. Lower right: 95MDF 95 percentile. B) Correlation of observed allele frequency (AF). 66/66 (100%) of expected variants were observed across both Manual and ALH-prepped samples. All samples were subsampled to 3M reads. (n=3 technical replicates per condition).

Conclusions

Our liquid reagent-based VARIANTPlex-HT workflow yielded libraries with parallel performance to our lyophilized products across a variety of input types, masses and panel sizes, including our most comprehensive panel to date, VARIANTPlex Complete Solid Tumor. The FUSIONPlex-HT workflow displayed greater numbers of unique reads supporting fusion calls across multiple input types, panels and at various dilutions of Fusion-positive material relative to lyophilized reagents. Automation of FUSIONPlex-HT and VARIANTPlex-HT workflows resulted in libraries with equivalent numbers of unique molecules and power to call variants relative to manually prepped samples. The FUSIONPlex-HT automated workflow detected all known fusions associated with FFPE input. Similarly, the VARIANTPlex-HT automated workflow detected all known variants associated with Myeloid input.

In conclusion, we have formulated and optimized automation-friendly liquid reagents that perform equivalently to legacy lyophilized products, providing a solution for mid- to high-throughput AMP-based NGS applications.

David Knupp, PhD
Research Scientist III (Integrated DNA Technologies)
dknupp@idtdna.com

