

ANTIBODY DISCOVERY

Synthetic gene and gene fragment solutions for drug development research



Reliable & fast turnaround times



High-quality dsDNA fragments



Easy to scale for automation

ANTIBODY RESEARCH IS CRUCIAL FOR ANTIBODY AND IMMUNO-ONCOLOGY DRUG DEVELOPMENT

The drug discovery process includes steps from target assessment to lead optimization and applies a variety of molecular engineering and synthetic biology methods. IDT has several gene and gene fragment solutions that can help with this early discovery phase.

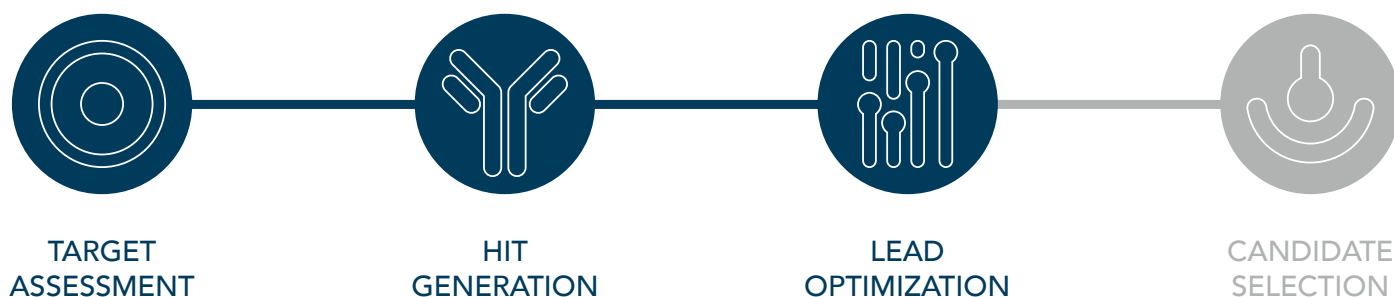


Figure 1. Antibody drug discovery workflow.

UTILIZING SYNTHETIC DNA FRAGMENTS FOR HIT GENERATION

Researchers can use several methods to generate antibodies. Synthetic DNA fragments are particularly useful when generating antibody fragments, such as single-domain nanobodies, which are easier to produce than traditional monoclonal antibodies, making them ideal for antibody-antigen interaction studies and to help identify potential hits for a desired drug target.

IDT offers three types of double-stranded DNA fragments: **gBlocks™ Gene Fragments**, **gBlocks HiFi Gene Fragments**, and **eBlocks™ Gene Fragments**. These fragments allow for convenient cloning of the antibody of choice and are also useful for a variety of *in vitro* selection techniques like phage or yeast display.

For Research Use Only. Not for use in diagnostic procedures.

> WWW.IDTDNA.COM

IDT GENE FRAGMENTS ARE A SCALABLE SOLUTION FOR HIGH-THROUGHPUT SCREENING

During the hit generation phase, hundreds of antibody variants may be generated and subsequently profiled. Identification of which hits have high titers and have reached pre-selection criteria for affinity, selectivity, and toxicity, is made easier with synthetic gene fragments from IDT (in particular, eBlocks Gene Fragments) as they are complementary to high-throughput screening approaches that may be used during this phase.

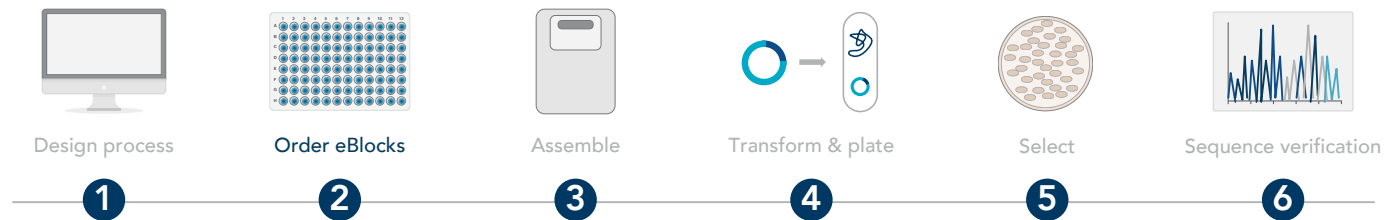


Figure 2. Automating high-throughput workflows. Researchers can design custom and ready-to-use gene fragments for DNA assembly, transformation, plating, and sequence verification steps for high-throughput screening in antibody discovery workflows.

SYNTHESIZING AND OPTIMIZING LEADS WITH DSDNA FRAGMENTS

Before selecting which leads to move forward as the therapeutic antibody candidate, promising leads can be further optimized. During the lead optimization phase, lead antibody molecules are selected, and additional alterations of the sequences can be done to improve efficacy and safety, as well as binding affinity. IDT **gene fragments** can be used to synthesize and optimize these sections (e.g., humanization via the Fc region or creation of bi-specifics), or site-directed mutagenesis libraries can be created. For researchers that prefer to skip cloning in-house, complete genes, which are cloned into vectors, can be **ordered**.

PRODUCT SUMMARY

| | eBlocks Gene Fragments | gBlocks Gene Fragments | gBlocks HiFi Gene Fragments | Gene Synthesis |
|---|------------------------|------------------------|-----------------------------|----------------|
| Key benefit | Speed | Quality | Premium quality | Ready-to-use |
| Available lengths | 300–1500 | 125–3000 | 1000–3000 | 25–5000+ |
| Median error rate | 1:5000 | 1:5000 | 1:12,000 | N/A* |
| Estimated shipping time (business days) | 1–3** | 2–8*** | 6–10 | 8+ |
| Yield | 200 ng | 250-1000 ng | 1000 ng | 4 µg |
| Format | Plate | Tube/Plate | Tube | N/A |

* Clonal genes contain no mutations present above IDT's sequencer-noise threshold

** Excluding Fridays as eBlocks plates that are resuspended and shipped on dry ice will not ship on Friday

*** This estimated shipping time is for tubes only. Plates estimated ship date is 10-15 business days (excluding Fridays)

RELATED PRODUCTS

| Antibody discovery workflow stage | IDT—related products |
|-----------------------------------|---|
| Target assessment | Alt-R™ Custom CRISPR gRNA Libraries, rhAmpSeq™ CRISPR Analysis System, xGen™ NGS products |
| Hit generation | oPools™ Oligo Pools, Ultramer™ DNA Oligonucleotides |
| Lead optimization | xGen™ NGS products, oPools Oligo Pools |

> FOR MORE INFORMATION, VISIT WWW.IDTDNA.COM/ANTIBODYDISCOVERY

For Research Use Only. Not for diagnostic procedures. Unless otherwise agreed to in writing, IDT does not intend these products to be used in clinical applications and does not warrant their fitness or suitability for any clinical diagnostic use. Purchaser is solely responsible for all decisions regarding the use of these products and any associated regulatory or legal obligations.

©2022 Integrated DNA Technologies, Inc. All rights reserved. gBlocks, eBlocks, Alt-R, rhAmpSeq, oPools, and Ultramer are trademarks of Integrated DNA Technologies, Inc. and registered in the USA. All other marks are the property of their respective owners. For specific trademark and licensing information, see www.idtdna.com/trademarks. Doc ID: RUO22-1151_001 10/22